



The HAPPYLEGS SYSTEM: Operation and Installation

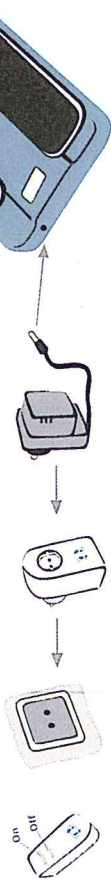

Happylegs improves the quality of life of all types of people who spend hours sitting or standing, almost immobile, and those who have reduced mobility.

- Extremely silent motor and low power consumption.
- Connects to the power supply with a transformer 7.5 V/9V/12V.*
- Three working speeds and remote control – On and Off.
- Height adjustable front legs.
- Highly durable product made in Spain with two years warranty.


* There are components for different voltages:
Europe 220 V, USA 110 V and UK 240 V.

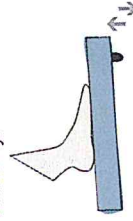

Ergonomic

Happylegs should be perfectly set on the floor. It has two front legs with adjustable height. It is advisable to place the apparatus with a slight inclination as this favours a better posture.


Always use the machine when seated. Use the remote control for turning the device "On" and "Off". Place your feet upon the two footpads. You should be relaxed and should let the machine move your feet without making any effort.

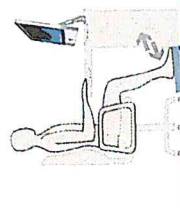

Select the required operational speed

It is important to check that the voltage is correct.

- 7.5 V: Slow Speed: Ideal for elderly people or those who use Happylegs for treatment of injuries.
- 9 V: Normal Speed: recommendable for continuous machine use. This should produce a pleasant feeling of well-being.
- 12 V: Fast Speed: recommendable when use of the machine is for periods of less than an two hours.

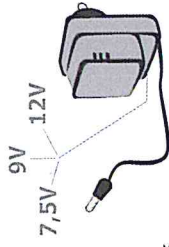

Cleaning and Maintenance


In order to use Happylegs correctly and for security reasons, the instructions outlined above should be followed.

We should never stand on top of the machine, as our weight would damage the internal operational mechanism. The machine should not be submitted to pressure for which it is not prepared.

Happylegs does not require any kind of maintenance and has been conceived to be used continuously for long periods of time. In order to clean it, a dry cloth should be used when the equipment is switched off.

Warning, see the ATTACHED DOCUMENTATION:

This machine is designed for use in all establishments, including homes and those directly connected to public low voltage power grid which supplies residential buildings. The quality of the power grid should be that of a typical commercial environment or hospital.

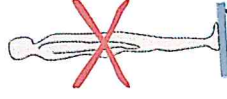

Transportation and Storage

This product is electrically operated; avoid spilling any liquid on Happylegs or storing in high humidity environments. Also avoid temperatures greater than 60° C, which could damage the equipment. Happylegs should be transported within its packaging. Avoid bumping or dropping and do not store weight greater than 40 k on it. If broken or damaged, contact Autoergon 2002, S.L. Not to be disassembled or repaired by unauthorised personnel.