## **Curera® Curved Pillow**

Pillow for torso support in seated position







## **Curera® Curved Pillow**

If a person makes too few spontaneous movements, there is an increased risk of pressure ulcers. Regular changes of position reduce the risk of pressure ulcers. Position as required in both reclining and seated position.

Curera® Positioning is a range of stable pillows, rolls, and wedges that meet most needs for support and relief. They are also a great help when changing dressings and for patient care.

Our soft hygiene covers and cotton covers are also available as accessories to the pillows and rolls.

The hygiene cover does not allow liquid or air to penetrate, but is vapour permeable (vapour is wicked away).

Make sure that the patient is positioned correctly so that airways are unobstructed.

Curera® Curved Pillow is curved for comfort, stability and safety in both the reclining position and seated position.



POSITIONING EXAMPLE

## **CLEANING**

The Hygiene Cover can be wiped clean with detergent and/or disinfectant. Machine wash at max. 95  $^{\circ}$ C.

Inner pillow: Machine wash at max. 95  $^{\circ}$ C. Can be tumble-dried or dried in an airing cupboard; avoid over-drying. For faster drying, run the spin cycle three or four times.

Cotton cover: Machine wash at max. 70 °C.


## **FACTS**

- Sizes: 85x65x30 cm
- Filling: 50 % HR foam/50 % polyester fibre
- Cover: Polyurethane/polyester.
- Accessories: Cotton cover.
- Warranty: 1 year
- Accessories:

Hygiene Cover: Polyurethane/polyester Cotton cover





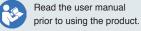

Care of Sweden and Curera are registered trademarks of Care of Sweden AB. © Care of Sweden AB, 2013. All rights reserved. Reproduction prohibited. Last updated: 09/04/2018