

**ECE Gruppe II - III**

15 - 36 kg

33.1 - 79.4 lb

# COMMANDER

The Specialist – Highly Adaptive Positioning  
Car Seat for children with special needs

  
Thomashilfen

# THOMASHILFEN COMMANDER

” The Thomashilfen Commander masters every challenge.

As a special seat for special requirements and positioning, it is ideally suited for highly individual supply solutions.

Children are supported wherever they need support and protection.

Especially for children with little upper body control, the fixed lateral trunk supports offer this extra support – the supports have proven themselves over many years of use in rehab strollers.

A huge selection of accessories and easy handling make the adjustment perfect. “



## Particularly suitable for:

- Cerebral Palsy (CP)
- Skull-brain trauma up to grade II
- Spina Bifida
- Rett Syndrome
- Autism
- Spinal Muscular Atrophy (SMA)

## Always recommended for


- Children with little upper body control



Thomashilfen has been providing children with special needs with rehab car seats for over 30 years. Today we also use our know-how to develop our own rehab car seats. All our many years of experience with the needs of children and parents flow into them.

A handwritten signature in black ink, reading "Gunnar Thomas".

Gunnar Thomas, owner

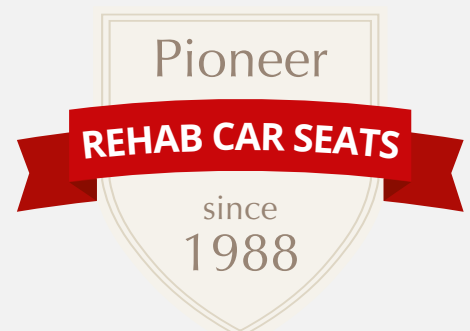

## REHAB-EQUIPMENT



\* optional

## UNIQUE COMBINATION



The unique, modular combination of **SWIVEL BASE** and **TILTING BASE** can be supplemented with a variety of accessories as required.



For easy handling, the release lever (draw latch) on the swivel base is easy to reach. It can be mounted either on the left or right front of the swivel base.

### UPHOLSTERY

cuddly upholstery for pleasant sitting comfort without annoying sweating, the breathable cover can be completely removed and washed at 40° C



### CUP HOLDER

removable large holder (right or left) to keep drinks, snacks and toys ready to hand



### SEAT TILT-IN-SPACE


with the integrated seat tilt-in-space +3° a resting position is quickly adjusted



## HEAD PADS

additional padding inserts can be added on the right and left side of the head pads, even better support for the child's head

  
**Thomashilfen**

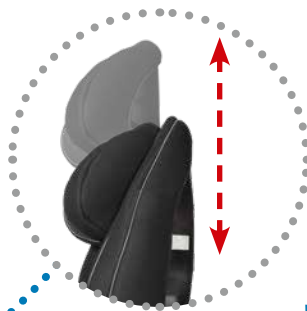

## HEAD SUPPORT

9-way height adjustable – easy and quick to adjust, grows step by step with your child



## SHOULDER BELT

"grows" with the child, without the need to rethread the belts



## BACK PLATE

sturdy, large back plate with hook-and-loop-fasteners for attaching the upper body supports – for optimum individual adjustment



## 5-POINT-BELT

integrated 5-point positioning belt with central adjuster for quick adjustment of the belt length, practical chest clip prevents twisting of the belts

Two versions **LATERAL TRUNK SUPPORTS** to choose from:



**FOAM CORE:** soft lateral trunk pads to attach directly onto the cover with hook-and-loop-fasteners



**METAL CORE:** stable lateral trunk supports for attaching to the back plate with hook-and-loop-fasteners

## ACCESSORIES – COMBINATION

Thomashilfen



**ADAPTER PLATE** –  
ALWAYS required for the  
following accessories



The unique modular combination  
of swivel base and tilting base can  
be combined with a variety of  
accessories.

**An adapter plate is ALWAYS  
required to use the accessories.**



**SWIVEL BASE** facilitates the  
transfer, easily accessible  
draw latch can be used left  
or right



**TILTING BASE** allows seat  
tilt in 3 steps (5.5°, 11°, 16°)  
for relaxed resting  
position



**TILTING BASE** and **SWIVEL  
BASE** combined for best  
possible handling and  
comfort



**FOOTREST WITH PAD**  
provides stable sitting posi-  
tion, children can firmly  
support their feet and do  
not slip out of the seat



**STAND, SEATFIX ADAPTER**  
and **STABILIZING BOW**  
guarantee stable seat  
support in the car



**SEAT DEPTH EXTENSION**  
allows the seat to grow  
with the child and ensures  
correct positioning of the  
thighs



**ABDUCTION BLOCK**  
(adjustable) ensures stable  
sitting position and pre-  
vents legs from crossing  
over

# THOMASHILFEN Commander




All details online:  
[www.thomashilfen.com/commander](http://www.thomashilfen.com/commander)



## Accessories



**ABDUCTION BLOCK** and **HIP SUPPORTS** support a stable sitting position

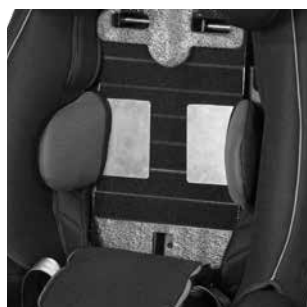

**LATERAL TRUNK SUPPORTS (METAL CORE)** can be positioned precisely on the back plate with hook-and-loop-fasteners for individual adjustment



**CROTCH PAD (LARGE)** protects the sensitive belly of the child, **SEAT WEDGE (BELOW)** allows +10° seat tilt



**SEAT WEDGE (INSIDE)** supports an upright pelvic position, **LATERAL TRUNK SUPPORTS (METAL CORE)** for optimum support of the upper body,

## Accessories



**TABLE** promotes the straightening up and stable posture of the child, provides the child with all-round support and safety



**ANTI-ESCAPE VEST** prevents the child from wriggling out of the belt



**CAR SEAT PROTECTION MAT WITH BACK** protects the seat upholstery and offers 3 practical compartments for toys at the front

A complete overview of accessories can be found on the order form on our website:  
[www.thomashilfen.com](http://www.thomashilfen.com)

Dimensions & weights	THOMASHILFEN Commander
<b>Seat depth</b>	34 cm / 13.4"
<b>Seat width</b> front (with supports)	33,5 cm (29 cm) / 13.2" (11.4")
<b>Seat width</b> back (with supports)	28,5 cm (27 cm) / 11.2" (10.6")
<b>Back height</b> (seat)	55 - 75 cm / 21.7 - 29.5"
<b>Shoulder height</b> (belt guide)	25 - 44 cm / 9.8 - 17.3"
<b>Seat tilt-in-space</b> (basic model)	3°
(with tilting plate)	5,5° / 11° / 16°
<b>Lower leg length</b>	23 - 40 cm / 9.1 - 15.7"
(Footrest directly attached to adapter plate)	
<b>Lower leg length</b>	28 - 45 cm / 11 - 17.7"
(Footrest attached to swivel base + tilting base)	
<b>Overall dimensions</b> (w x h x d)	55 x 65 x 45 cm / 21.7 x 25.6 x 17.7"
<b>Weight of seat</b>	approx. 8,5 kg / 18.7 lb
<b>max. user weight</b>	15 - 36 kg / 33.1 - 79.4 lb



5 upholstery designs available:



### Basic configuration:

Seat and back with integrated 9-fold height-adjustable headrest; back plate (with hook-and-loop fasteners); seat cushion; back cushion; head cushion; additional head cushion inserts; crotch pad, small; 5-point positioning belt with central adjuster; chest clip; shoulder protector; cup holder (can be used right or left); removable covers, machine washable up to 40° C



Thomas Hilfen für Körperbehinderte GmbH & Co. Medico KG · Walkmühlenstraße 1 · D - 27432 Bremervörde  
 Phone: +49 (0)4761 886-63 or -68 · Fax: +49 (0)4761 886-19 · [www.thomashilfen.com](http://www.thomashilfen.com)