

'RE FOR TESTING AND CERTIFICATION - MECH-TEST

Mechanical Laboratory

05-077 Warszawa-Wesoła, ul. Klonowa 22 tel.: +48 603 23-26-45, e-mail: cbc.mech.test@gmail.com, www.cbc.org.pl

Date 5.01.2016

TEST REPORT NO. *CBC* –148/2015

Page 1 of 4

Subject of testing:

Walking aids with built-in handgrips and three or

Classification according to PN-EN ISO 9999:2011:12 06 06

more legs of which two or more are having wheels,

which provide support whilst walking

Type / Model:

Explorer Rollator

LARGE

Nr kat.: ROLOOE

Number of specimens: 1

SN: 0001

Manufacturer:

REHASENSE Sp. z o.o.

Suleiowska 45

97-300 Piotrków Trybunalski

Applicant:

A-Net s.c.

93-469 Łódź.

ul. Łaskowice174

Kind of testing

Mechanical testing for conformity with PN-EN ISO 11199-2 : 2005

excluding clause 4.10, 6.2, 6.3

Test started: 17.12.2015

Test finished: 5.01.2016

Approved by:

inż. Andrzej Tkaczyk

Special comments / enclosures:

Copyright © 2012 by Centre for Testing and Certification (applicable to report form)

Test results refer only to tested units.

Test results reported here are not applicable to the further modifications of the product affecting its structure, material or technology.

This test report shall be neither copied differently as in the whole nor be published without written consent of the Laboratory.



Page: 2 of 4

CHARACTERISTIC OF PRODUCT

Name: LARGE


Dimension of rollator: -
Nr. kat: ROL00E

Product code: --

Maximum permissible user mass:200 kgMass of rollator:9,365 kgDescriptionPHOTO OF PRODUCT

Description					
Ele	ments/parameters/materials/di	mensions	Co	omments	
Dimensions od walking rollator (fig. 2 PN-EN ISO 11199-2)	Distance between	560 mm			
	handgrips (dimension 2)				
	Angle between of handgrip	00			
	axis and direction of				
	movement (α)				
	Height of rollator	907 mm	mi	in.	
	(dimension 6)	1121 mm	max.		
P &	Width of rollator	680 mm			
ns .	(dimension 5)				
sio	Turning width	950 mm			
Ten (fig	(dimension 1)				
iii	Length of rollator	760 mm	T		
	(dimension 4)				
Dimen	nsions of folded rollator (mm) 760 x 909		25.	2	
	Handgrip - diameter	37 mm		Anatomical	
Fig. 3				handgrip	
	Handgrip - length	120 mm	******		
	Front wheels- quantity	2		castor	
ato.	Front wheels - diameter	282 mm		wheels	
olla	Front wheels – width	43 mm			
f r	Front wheels - brake	none			
Wheels of rollator	Rear wheels - quantity	2			
	Rear wheels - diameter	232 mm			
₩	Rear wheels - width	44 mm			
	Rear wheels - brake	Included			
Tip	Diameter				
	Material	Not any			
	Colour				
Material of rolla- tor (fig. 1)	Front legs	Aluminum,			
	Bracing member (no. 8)	steel, plastic,			
	Rear legs	Bolts, nuts			
100	Height adjusting device (no.				
erria fig.	4)				
Material o tor (fig. 1)	Handgrip (no 5),	plastic			
 	Brake elements				



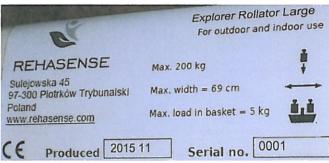











Report no.: CBC-148/2015

Page: 3 of 4

men ding se	0	1						
Requiremen ts according to clause	Test me- thod according to clause		charac	Checked characteristics/assemblies/parameters		Real	Test result	Comments
4.1	Measur.	Manoeuvrability				ø 282 mm width 41 mm Conf.	Pos.	ø front wheels ≥75mm <u>outdoor intended rollator:</u> ø front wheels ≥180mm width of wheels ≥22mm
4.2	5.3	Forward-direction stability				19,5 ° Conf.	Pos.	Stability required ≥ 15°
4.2	5.4	Backward-direction stability				10,5 ° Conf.	Pos.	Stability required ≥ 7°
4.2	5.5	Sideway-direction stability Stability – forwards			4,1 ° Conf.	Pos.	Stability required ≥ 3,5°	
4.2	5.6			akat haa dein ayyyaan aylinday	forwards	18,0° Conf.	Pos.	Stability required ≥ 15°
		with loaded basket, bag, drip, oxygen cylinder backwards side		14,8 ° Conf. 5,5 ° Conf.	Pos. Pos.	Stability required $\geq 7^{\circ}$ Stability required $\geq 3.5^{\circ}$		
4.3	V/I		Servicing wheels	facility during rollator motion w		Conf.	Pos.	Stability required 23,5
	V/I		resting se	rakes in rollator with more than 2 at or intended for outdoor use	2 wheels and	Conf.	Pos.	
	5.7.1.1			o distance (fig. 4, dimension 1)		65 mm Conf.	Pos.	≤ 75 mm
-	5.7.1	F F		brake effectiveness		Conf.	Pos.	Movement of rollator ≤10 mm in 1 minute
L	Measur.	- L		et parking brake		30N Conf.	Pos.	≤ 60 N
	Measur.	1	Force to r	elease parking brake		10N Conf.	Pos.	≤ 40 N
	5.7.2			rake effectiveness		Conf.	Pos.	Movement of rollator ≤ 10 mm in 1 minute
	V/I		Possibility	y to compensate brake wear		Conf.	Pos.	
	V/I			adversely affected by folding, unactions of rollator	nforlding or	Conf.	Pos.	
4.4	Measur. V/I	Handgrip				37 mm Conf.	Pos.	Width of handgrip ≥20mm and ≤50mm
4.5	Measur. V/I	Leg section and tip				N/A	ø tip ≥35mm (tested rollator is equipped in four wheels)	
4.6	5.10	Resting seat – static loading durability				Conf.	Pos.	1 minute under load 1,2 x user's weight±2% (240kg)
4.7	5.12	Mechanical durability		Fatigue test		Conf.	Pos.	200 000 cycles with load. 160kg±2%, f=1Hz
4.7	5.11			Static loading test	6. de	Conf.	Pos.	loading 240kg±2%, 5sek.
4.8	V/I		Adjusting devices			Conf.	Pos.	
4.9	5.14 ISO 10993-	Foldi	ng mecha	nism ocompatibility of material with h	uman hadai	Conf.	Pos.	
4.11	1 V/I	ials	DI				N/T	
-		Materials and finish	wi	ee of discolouring of skin or cloth th rollator materials	ing in contact	Conf.	Pos.	
	V/I	-	Bu	errs, shar edges, projections Marking and labe	lling of product	Conf.	Pos.	
6.2	V/I	a) Ma	aximum u		ining of product		N/T	T
				afe working load (SWL) to be ma	rked on accessori		N/T	
		 c) Maximum allowed angle between the longitudinal centreline of handle and the direction of motion, if the handles are sideways adjustable 				N/T		
	ļ	d) Ma	anufacture	r's name or trade name and addre			N/T	
	[r's model identification name and	d/or number		N/T	
		g) Ma	aximum ex	ear of manufacture xtension of the height adjustment	, marked on the		N/T N/T	
	+		adjusting members n) Maximum width of the rollator					
	}			ided for outdoor/indoor use			N/T	
4.10	V/I	Warning showing allowed angle between handle axis and direction movement or physical stop of angle adjusting			1 of	N/T N/T		



Mechanical Laboratory of CBC			Report no.: CBC-148/2015 Page : 4 of 4				
Contents of user manual and/or assembly manual or clear and indelible marking of product							
6.3	V/I	a) Maximum rollator height		N/T			
1		b) Minimum rollator height		N/T			
		c) maintenance and cleaning instructions, including a description of the method and suitable cleaning agents and any precautions needed to avoid corrosion and/or ageing of the materials used in construction of the rollator		N/T			
		d) Instructions for assembly, adjustment of all kinds, folding and unfolding		N/T			
		e) Warnings and advice about precautions relating to safe distances between moving and stationary parts (see EN 12182, Clauses 12 and 13, for guidance)	_	N/T			
		f) Maximum safe working load (SWL) for load carrying accessories such as basket, tray, shopping bag, etc.		N/T			
4.10	V/I	Warning in user manual on consequences of such an adjustment of angle between handle longitudinal axis and direction of movement outside allowed value (when handles are adjustable aside).		N/T			
		TEST CONDITIONS					
Ambient temperature		19°C		Required temperature 21°C ±5°C			
Relative humidity of air:			55 %		Not required		
Commer	nts:						
		maximum height adjustment of rollator.					
		he least stabble position of self-adjusting wheels.					
Tests per	formed with ha	ndles positioned at their maximum (allowed) angle to the direction of moti	on (when a	djustment is	s possible).		
		ity test, static loading test, fatigue test.					
	tor was tested.						
		on before testing any visible defects that could have influence on te	est results	were not st	tated.		

Pos. – positive; Neg – negative; N/T – not tested; N/A – not applicable; N/R – not required , N/O – not occurred , V/I. – visual inspection, Conf. – conformed.

NOTE 1: Deformation – 58 mm, elastic deformation – 54 mm, permanent deformation – 4 mm.

CONCLUSIONS:

Testing object **conforms** with requirements of PN-EN ISO 11199-2: 2005, excluding clause 4.10, 6.2, 6.3 in scope of mechanical testing ordered by client, excluding biocompatibility tests of material with human body according to PN-EN ISO 10993-1:2010.

---- END ----

